

# One Plug. One Size.

## SNUGPLUG™


### That's *all* you need!

SNUGPLUG™ is the result of our extensive research on Dry Eye syndrome over the last 30 years. It relieves the discomfort associated with Dry Eye symptoms quickly and effectively to improve the patient's quality of life.

- One-size-fits-most punctum plug **eliminates sizing** and **simplifies stocking**
- Pre-stretched on inserter **avoids the need for dilation** and **facilitates insertion**
- Returns to its natural shape when released from inserter
- 1.1 mm bulb assures snug fit and **eliminates pop-out**
- Soft collarette **provides patient comfort**




Visit [FCI-Ophthalmics.com](http://FCI-Ophthalmics.com) to watch the SnugPlug™ informational video.

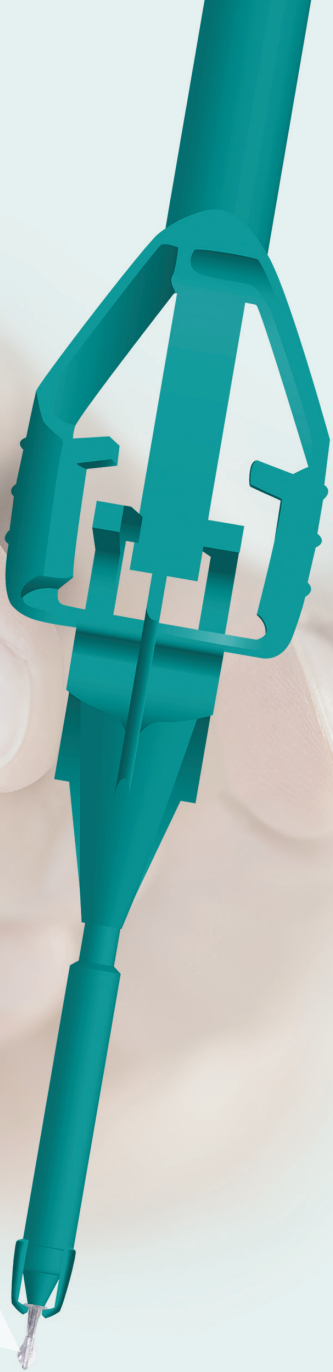

## PERMANENT PUNCTAL PLUGS

### SNUGPLUG™

- Indicated for punctal occlusion
- One-size-fits-most
- No sizing needed
- High retention rate\*
- No dilation needed
- Ultra-thin collarette
- Medical grade silicone
- Individually packaged, sterile

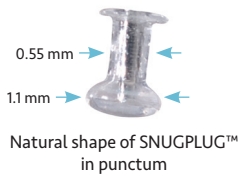

SNUGPLUG™	SIZE	PACKAGING
S2.4002u	One-size-fits-most	Box of 2

All FCI products are completely LATEX-FREE.

FCI products are sold worldwide through a network of independent distributors and representatives. To place an order or find your distributor or representative, call toll-free 800.932.4202.

#### Payment Terms

30 days net  
Shipping prepaid and added to invoice

#### Returns

Defective merchandise may be returned for replacement only. Merchandise shipped in error by FCI/FCI Distributor may be returned for credit or replacement within 30 days from date of invoice. Note: Some products are non-returnable.

- Prior authorization **MUST** be obtained.
- Returned merchandise must be properly packaged to prevent damage during shipping.
- Returned merchandise must be in original condition – unopened, unmarked and undamaged.

Restocking charge may apply.

For full return policy, visit [fci-ophthalmics.com/company/terms-and-conditions](http://fci-ophthalmics.com/company/terms-and-conditions)

#### Guarantee

FCI products are unconditionally guaranteed against manufacturing defects. Any product found to have a manufacturing defect will be replaced free of charge.

The products listed in this publication may not be registered in your country. Contact your local distributor to find out what is available to you.

This publication is intended for distribution in the USA and Canada only.

## CONTACT FCI

#### Orders

[orders@fci-ophthalmics.com](mailto:orders@fci-ophthalmics.com)

#### Telephone

800.932.4202  
781.826.9060

#### Fax

781.826.9062

#### Information

[info@fci-ophthalmics.com](mailto:info@fci-ophthalmics.com)



## JOIN OUR OPHTHALMIC COMMUNITY



LINKEDIN



YOUTUBE



EYETUBE



# Dry Eye Catalog

**FCI-Ophthalmics.com**  
800.932.4202

**FCI-Ophthalmics.com**



FCI Ophthalmics  
30 Corporate Park Drive, Suite 310/320  
Pembroke, MA 02359

For more information, please contact your local representative:



\* Javate RM, Dy IE, Buyucan KF, Ma Guerrero EE. Retention rates and benefits of Painless Punctal Plug F(TM) in dry eye patients. Orbit. 2016 Jun;35(3):126-31. doi: 10.1080/01676830.2016.1176204. Epub 2016 May 10. PMID: 27163906.


PERMANENT PUNCTAL PLUGS

PRELOADED PUNCTAL PLUGS®

FCI's PRELOADED PUNCTAL PLUGS® are "ready-set plugs" with an ultra-thin silicone collarette for maximum patient comfort. For ease of use, FCI punctal plugs are preloaded on a single use inserter/dilator and are easily removed with forceps.

Designed by J.A. Bernard, M.D.

- Indicated for punctal occlusion
- Original design
- Ergonomic slanted collarette
- Preloaded on single use inserter/dilator
- Good tolerance and optimal patient comfort
- Medical grade silicone
- Individually packaged, sterile

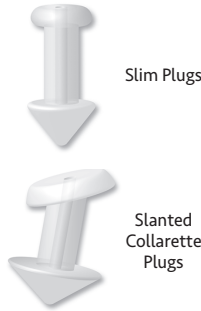


**Slim Plugs**  
The Slim Plugs offer a flat collarette and a slightly smaller bulb for easier insertion into small puncta.

**Exclusive Slanted Collarette Plugs**  
The Extra Small to Large plugs feature the unique slanted collarette that hugs the natural curve of the lid margin. The slightly larger bulb resists migration and is designed to prevent pop-out.

SLIM PLUGS FLAT COLLARETTE	SIZE	WIDTH	PACKAGING
<b>S2.3632</b>	Slim Mini	0.4 mm	Box of 2
<b>S2.3642</b>	Slim Petite	0.5 mm	Box of 2

SLANTED COLLARETTE PLUGS	SIZE	WIDTH	PACKAGING
<b>S2.3652</b>	Extra Small	0.6 mm	Box of 2
<b>S2.3672u</b>	Mini	0.7 mm	Box of 2
<b>S2.3422</b>	Small	0.8 mm	Box of 2
<b>S2.3122</b>	Regular	0.9 mm	Box of 2
<b>S2.3522</b>	Large	1.0 mm	Box of 2

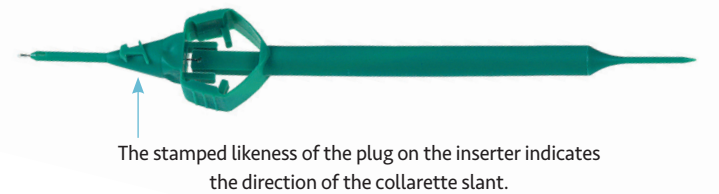

TEARSAVER PLUGS

The tapered shaft design helps to anchor the plug securely in the punctum for less risk of pop-out.

Preloaded on a single use inserter/dilator for ease of insertion.

- Indicated for punctal occlusion
- Medical grade silicone
- Permanent yet removable
- Individually packaged, sterile



TEARSAVER PLUGS	SIZE	PACKAGING
<b>814004</b>	0.4 mm	Box of 2
<b>814005</b>	0.5 mm	Box of 2
<b>814006</b>	0.6 mm	Box of 2
<b>814007</b>	0.7 mm	Box of 2
<b>814008</b>	0.8 mm	Box of 2

PVP PERFORATED PLUGS

The plugs are coated with a thin layer of polyvinylpyrrolidone (PVP)\* to improve the hydrophilic properties of silicone. This enhances the efficiency of tear drainage and limits the accumulation of debris at the surface.

Designed by J.A. Bernard, M.D.

- Indicated for partial occlusion or punctal stenosis
- PVP coated medical grade silicone
- Pre-loaded on a single use plug inserter to facilitate insertion
- Slanted collarette for a better anatomical fit against the eyelid margin
- Individually packaged, sterile



PERFORATED PLUGS	SIZE	WIDTH	CHANNEL DIAMETER	PACKAGING
<b>S1.3522</b>	Mini	0.7 mm	0.4 mm	Box of 2
<b>S1.3512</b>	Medium	0.9 mm	0.6 mm	Box of 2



Perforation through length of plug acts as a conformer, allowing tears to drain while maintaining punctal opening.

\* Th Malet, B Challier, N David, A Bertrand, J-L George "Clinical and scintigraphic comparison of silicone and polyvinylpyrrolidone coated silicone PERFORATED PLUGS". Br J Ophthalmol 1998;82:1416-1419.

NON-STERILE PLUGS

FCI's NON-STERILE PUNCTAL PLUGS® are "ready-set plugs" pre-mounted on a wheel for quick and easy loading onto a plug inserter.

- Indicated for punctal occlusion
- Medical grade silicone
- Includes 2 Single Use Punctum Dilator and Plug Inserters (\$1.3090)
- Non-sterile

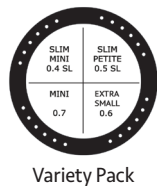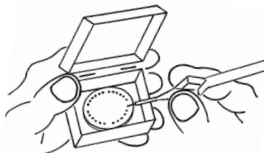

SLIM PLUGS FLAT COLLARETTE	SIZE	WIDTH	PACKAGING
<b>S2.3324u</b>	Slim Mini	0.4 mm	Box of 24
<b>S2.3424u</b>	Slim Petite	0.5 mm	Box of 24

SLANTED COLLARETTE PLUGS	SIZE	WIDTH	PACKAGING
<b>S2.3524u</b>	Extra Small	0.6 mm	Box of 24
<b>S2.3724u</b>	Mini	0.7 mm	Box of 24

VARIETY PACK	SIZE	WIDTH	PACKAGING
<b>S2.3624u</b>	6 plugs of each size	0.4, 0.5, 0.6, 0.7 mm	Box of 24



DRY EYE ACCESSORIES

REUSABLE DOUBLE DILATOR

- Stainless steel
- Sterile

REUSABLE DOUBLE DILATOR	TAPER	PACKAGING
<b>A8-4050</b>	Fine & Medium	Box of 1



TEMPORARY PUNCTAL PLUGS

COLLAGEN PLUGS

- Indicated for use in the diagnosis of Dry Eye
- To assess the effectiveness of permanent punctal occlusion
- Dissolves within 7 to 10 days
- 2 mm length
- Sterile

6 plugs per packet, 12 packets per box (72 total)

Variety Pack: 6 plugs per packet, 4 packets of each size per box (72 total)



COLLAGEN PLUGS	SIZE	PACKAGING
<b>CP-0002</b>	0.2 mm	12 Packets of 6
<b>CP-0003</b>	0.3 mm	12 Packets of 6
<b>CP-0004</b>	0.4 mm	12 Packets of 6
<b>CP-0234</b>	Variety Pack	12 Packets of 6

PROLONG™ ABSORBABLE LONG TERM PLUGS

- Indicated for the treatment of Dry Eye syndrome and Dry Eye conditions due to various ocular surface diseases
- For post-surgical Dry Eye symptoms
- Lasts up to 3 months
- 2 mm length
- Sterile

2 plugs per packet, 10 packets per box (20 total)

Variety Pack: 3 packets 0.3 mm, 5 packets 0.4 mm, 2 packets 0.5 mm (20 total)



PROLONG™ PLUGS	SIZE	PACKAGING
<b>TP-2000</b>	0.2 mm	10 Packets of 2
<b>TP-3000</b>	0.3 mm	10 Packets of 2
<b>TP-4000</b>	0.4 mm	10 Packets of 2
<b>TP-5000</b>	0.5 mm	10 Packets of 2
<b>TP-3450</b>	Variety Pack	10 Packets of 2

ULTRAPLUG™ EXTENDED WEAR ABSORBABLE PLUGS

- Indicated for use in the treatment of Dry Eye syndrome as a result of seasonal allergies, eye-related surgery, or contact lens wear
- E-Caprolactone-L-Lactide copolymer (PCL) – essentially absorbed within 2-6 months
- 2 mm length
- Sterile

2 plugs per packet, 10 packets per box (20 total)



ULTRAPLUG™	SIZE	PACKAGING
<b>EPLUG2</b>	0.2 mm	10 Packets of 2
<b>EPLUG3</b>	0.3 mm	10 Packets of 2
<b>EPLUG4</b>	0.4 mm	10 Packets of 2
<b>EPLUG5</b>	0.5 mm	10 Packets of 2