

SRS II

Self-Retaining Bicanaliculus Intubation Set II

INDICATIONS

For Horizontal Lacrimal
Duct Obstruction and
Punctal Stenosis



FCI®

a ZEISS company



www.FCI-Ophthalmics.com

SELF-RETAINING BICANALICULUS INTUBATION SET II (SRS II)

The SRS II enables intubation of the upper and lower canaliculi without having to enter the nasal passageway. The 0.64 mm wide tube is made of medical grade silicone and is available in three lengths – 25 mm, 30 mm and 35 mm. Each end of the tube has an anchor-shaped head with two flexible winglets that fold inward during insertion through the punctum. After passage through the common canaliculus, the winglets open in the lacrimal sac for secure fixation. A centrally placed marking on the tube acts as a reference point and allows verification of proper stent positioning.

FEATURES & BENEFITS

QUICK & SIMPLE


- In-office procedure using anesthetic eye drops
- Less invasive
- Virtually painless for the patient

EASY & SECURE

- Flexible winglet design enables secure fixation
- Guides give rigidity to the tube and facilitate placement of the stent
- Central mark offers control of proper stent positioning

CONVENIENT

- Available in three lengths – 25 mm, 30 mm and 35 mm
- Includes disposable dilator
- Disposable Sizer (S1-1289) sold separately

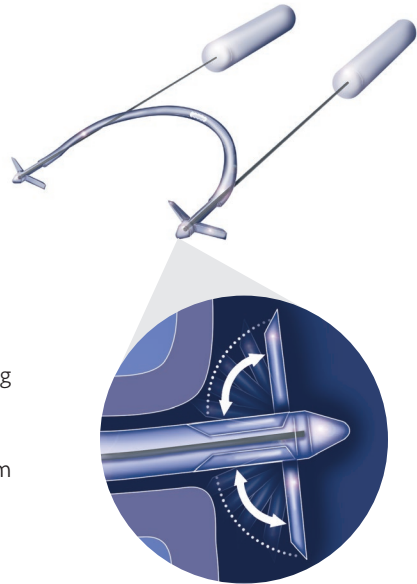

Insertion Through Punctum:
Winglets Fold Inward

Passage Through Canaliculus:
Winglets Open in Lacrimal Sac

S1.1390	25 mm	Box of 1, sterile
S1.1391	30 mm	Box of 1, sterile
S1.1392	35 mm	Box of 1, sterile
S1.1289	Disposable Sizers	Box of 5, sterile



www.FCI-Ophthalmics.com



FCI OPHTHALMICS

30 Corporate Park Drive, Suite 310/320
Pembroke, MA 02359

EMAIL

info@fci-ophthalmics.com

TELEPHONE

800.932.4202 • 781.826.9060

FAX

781.826.9062