



# Dry Eye Catalog

[FCI-Ophthalmics.com](http://FCI-Ophthalmics.com)



**FCI**®

a ZEISS company


# One Plug. One Size.

## SNUGPLUG®

### That's *all* you need!

SNUGPLUG® is the result of our extensive research on Dry Eye syndrome over the last 30 years. It relieves the discomfort associated with Dry Eye symptoms quickly and effectively to improve the patient's quality of life.

- One-size-fits-most punctum plug **eliminates sizing** and **simplifies stocking**
- Pre-stretched on inserter **avoids the need for dilation** and **facilitates insertion**
- Returns to its natural shape when released from inserter
- 1.1 mm bulb assures snug fit and **eliminates pop-out**
- Soft collarete **provides patient comfort**

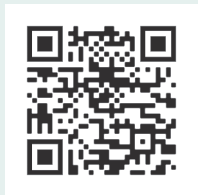


Visit [FCI-Ophthalmics.com](http://FCI-Ophthalmics.com) to watch the SnugPlug® informational video.

## PERMANENT PUNCTAL PLUGS

### SNUGPLUG®

- Indicated for punctal occlusion
- One-size-fits-most
- No sizing needed
- High retention rate\*
- No dilation needed
- Ultra-thin collarete
- Medical grade silicone
- Individually packaged, sterile

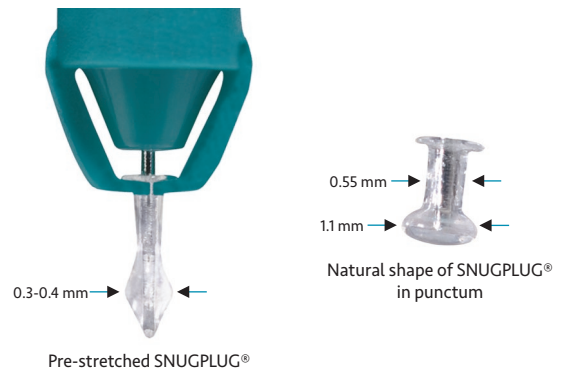

SNUGPLUG®	SIZE	PACKAGING
<b>S2.4002u</b>	One-size-fits-most	Box of 2

\* Javate RM, Dy IE, Buyucan KF, Ma Guerrero EE. Retention rates and benefits of Painless Punctal Plug F(TM) in dry eye patients. Orbit. 2016 Jun;35(3):126-31. doi: 10.1080/01676830.2016.1176204. Epub 2016 May 10. PMID: 27163906.


# PERMANENT PUNCTAL PLUGS

## PRELOADED PUNCTAL PLUGS

FCI's PRELOADED PUNCTAL PLUGS are "ready-set plugs" with an ultra-thin silicone collarette for maximum patient comfort. For ease of use, FCI punctal plugs are preloaded on a single use inserter/dilator and are easily removed with forceps.

**Designed by J.A. Bernard, M.D.**

- Indicated for punctal occlusion
- Original design
- Ergonomic slanted collarette
- Preloaded on single use inserter/dilator
- Good tolerance and optimal patient comfort
- Medical grade silicone
- Individually packaged, sterile

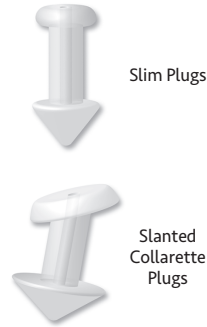

### Slim Plugs

The Slim Plugs offer a flat collarette and a slightly smaller bulb for easier insertion into small puncta.

### Exclusive Slanted Collarette Plugs

The Extra Small to Large plugs feature the unique slanted collarette that hugs the natural curve of the lid margin. The slightly larger bulb resists migration and is designed to prevent pop-out.

SLIM PLUGS FLAT COLLARETTE	SIZE	WIDTH	PACKAGING
<b>S2.3632</b>	Slim Mini	0.4 mm	Box of 2
<b>S2.3642</b>	Slim Petite	0.5 mm	Box of 2

SLANTED COLLARETTE PLUGS	SIZE	WIDTH	PACKAGING
<b>S2.3652</b>	Extra Small	0.6 mm	Box of 2
<b>S2.3672u</b>	Mini	0.7 mm	Box of 2
<b>S2.3422</b>	Small	0.8 mm	Box of 2
<b>S2.3122</b>	Regular	0.9 mm	Box of 2
<b>S2.3522</b>	Large	1.0 mm	Box of 2



The stamped likeness of the plug on the inserter indicates the direction of the collarette slant.






## TEARSAVER® PLUGS

The tapered shaft design helps to anchor the plug securely in the punctum for less risk of pop-out.

- Preloaded on a single use inserter/dilator for ease of insertion
- Indicated for punctal occlusion
- Medical grade silicone
- Permanent yet removable
- Individually packaged, sterile



TEARSAVER® PLUGS	SIZE	PACKAGING
<b>S2.3842</b>	0.4 mm	Box of 2
<b>S2.3852</b>	0.5 mm	Box of 2
<b>S2.3862</b>	0.6 mm	Box of 2
<b>S2.3872</b>	0.7 mm	Box of 2
<b>S2.3882</b>	0.8 mm	Box of 2

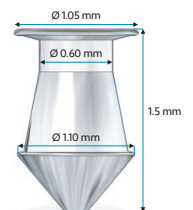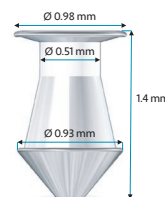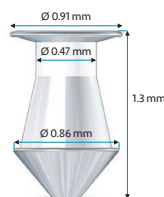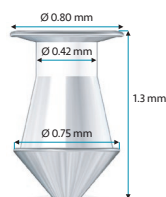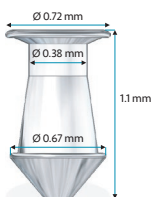

## PVP PERFORATED PLUGS

The plugs are coated with a thin layer of polyvinylpyrrolidone (PVP)\* to improve the hydrophilic properties of silicone. This enhances the efficiency of tear drainage and limits the accumulation of debris at the surface.

Designed by J.A. Bernard, M.D.

- Indicated for partial occlusion or punctal stenosis
- PVP coated medical grade silicone
- Preloaded on a single use plug inserter to facilitate insertion
- Slanted collarette for a better anatomical fit against the eyelid margin
- Individually packaged, sterile



PERFORATED PLUGS	SIZE	WIDTH	CHANNEL DIAMETER	PACKAGING
<b>S1.3522</b>	Mini	0.7 mm	0.4 mm	Box of 2
<b>S1.3512</b>	Medium	0.9 mm	0.6 mm	Box of 2



Perforation through length of plug acts as a conformer, allowing tears to drain while maintaining punctal opening.

\* Th Malet, B Challier, N David, A Bertrand, J-L George "Clinical and scintigraphic comparison of silicone and polyvinylpyrrolidone coated silicone PERFORATED PLUGS". Br J Ophthalmol 1998;82:1416-1419.

## NON-STERILE PLUGS

FCI's NON-STERILE PUNCTAL PLUGS are pre-mounted on a wheel for quick and easy loading onto a plug inserter.

- Indicated for punctal occlusion
- Medical grade silicone
- Includes 2 single use Punctum Dilator and Plug Inserters (S1.3090)
- Non-sterile

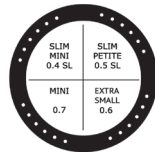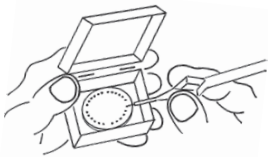

Variety Pack

SLIM PLUGS FLAT COLLARETTE	SIZE	WIDTH	PACKAGING
<b>S2.3324u</b>	Slim Mini	0.4 mm	Box of 24
<b>S2.3424u</b>	Slim Petite	0.5 mm	Box of 24

SLANTED COLLARETTE PLUGS	SIZE	WIDTH	PACKAGING
<b>S2.3524u</b>	Extra Small	0.6 mm	Box of 24
<b>S2.3724u</b>	Mini	0.7 mm	Box of 24

VARIETY PACK	SIZE	WIDTH	PACKAGING
<b>S2.3624u</b>	6 plugs of each size	0.4, 0.5, 0.6, 0.7 mm	Box of 24



## DRY EYE ACCESSORIES

### REUSABLE DOUBLE DILATOR

- Stainless steel
- Non-sterile

REUSABLE DOUBLE DILATOR	TAPER	PACKAGING
<b>A8-4050</b>	Fine & Medium	Box of 1



## TEMPORARY PUNCTAL PLUGS

### FCI PROLONG® Dx

- Indicated for use in the diagnosis of Dry Eye
- To assess the effectiveness of permanent punctal occlusion
- Dissolves within five days
- 2 mm length
- Sterile
- Made of collagen



6 plugs per packet, 10 packets per box (60 total)

Variety Pack: 6 plugs per packet – 3 packets of 0.2 mm, 3 packets 0.3 mm, 4 packets of 0.4 mm (60 total)



FCI PROLONG® Dx	WIDTH	PACKAGING
<b>S2.5202</b>	0.2 mm	10 Packets of 6
<b>S2.5203</b>	0.3 mm	10 Packets of 6
<b>S2.5204</b>	0.4 mm	10 Packets of 6
<b>S2.5224</b>	Variety Pack	10 Packets of 6

### FCI PROLONG®

- Indicated for the treatment of Dry Eye syndrome and Dry Eye conditions due to various ocular surface diseases
- For post-surgical Dry Eye symptoms
- Absorbs within 90 days
- 2 mm length
- Sterile



2 plugs per packet, 10 packets per box (20 total)

Variety Pack: 2 plugs per packet – 3 packets 0.3 mm, 5 packets 0.4 mm, 2 packets 0.5 mm (20 total)



FCI PROLONG®	WIDTH	PACKAGING
<b>S2.5002</b>	0.2 mm	10 Packets of 2
<b>S2.5003</b>	0.3 mm	10 Packets of 2
<b>S2.5004</b>	0.4 mm	10 Packets of 2
<b>S2.5005</b>	0.5 mm	10 Packets of 2
<b>S2.5035</b>	Variety Pack	10 Packets of 2

### FCI PROLONG® +

- Indicated for use in the treatment of Dry Eye syndrome as a result of seasonal allergies, eye-related surgery, or contact lens wear
- E-Caprolactone-L-Lactide copolymer (PCL)
- Absorbs within 180 days
- 2 mm length
- Sterile



2 plugs per packet, 10 packets per box (20 total)

Variety Pack: 2 plugs per packet – 3 packets 0.3 mm, 4 packets 0.4 mm, 3 packets 0.5 mm (20 total)



FCI PROLONG® +	WIDTH	PACKAGING
<b>S2.5103</b>	0.3 mm	10 Packets of 2
<b>S2.5104</b>	0.4 mm	10 Packets of 2
<b>S2.5105</b>	0.5 mm	10 Packets of 2
<b>S2.5135</b>	Variety Pack	10 Packets of 2

## CONTACT FCI

### Orders

orders@fci-ophthalmics.com

### Telephone

800.932.4202

781.826.9060

### Fax

781.826.9062

### Information

info@fci-ophthalmics.com



## JOIN OUR OPHTHALMIC COMMUNITY



LINKEDIN



YOUTUBE



EYETUBE

This publication is intended for distribution in North America only.

# FCI-Ophthalmics.com



For more information, please contact your local representative: